
PL

EN

CS

DE

NL
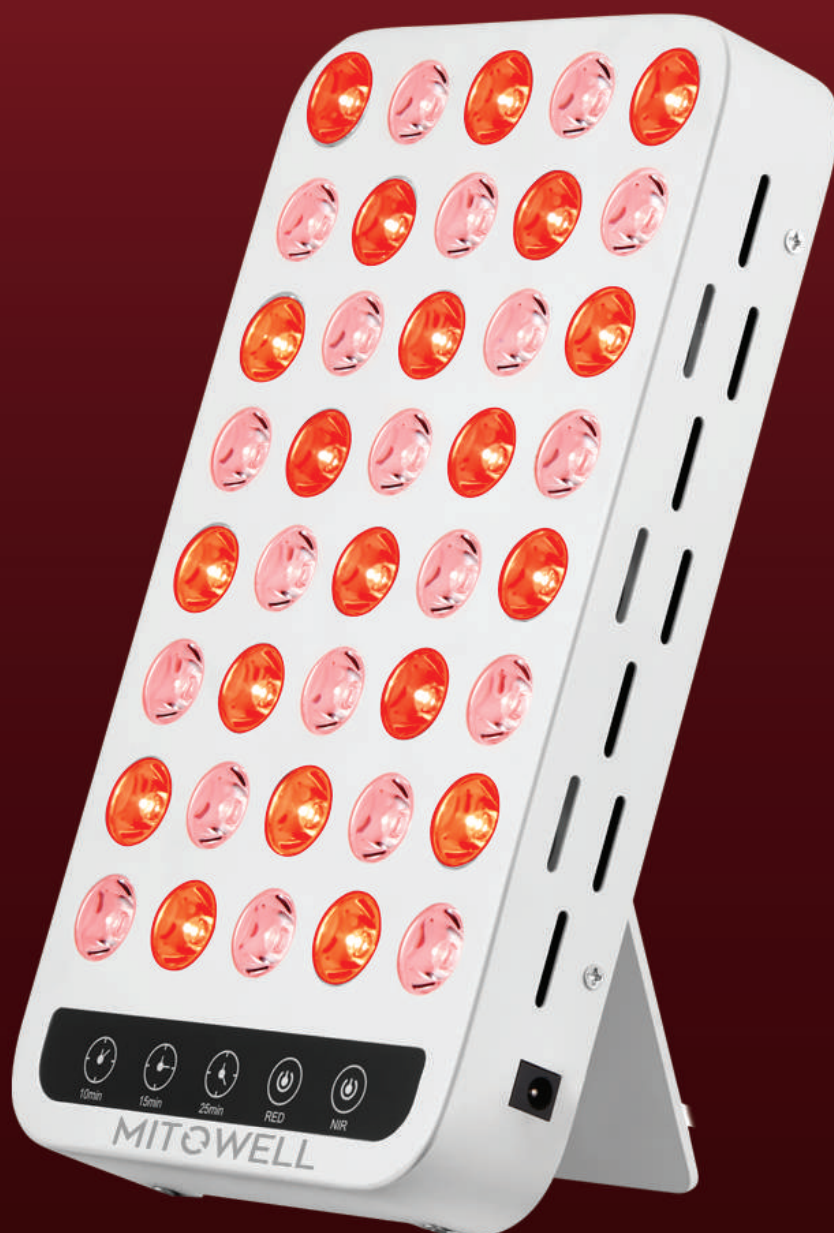
# MITOWELL

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M200

# SPIS TREŚCI

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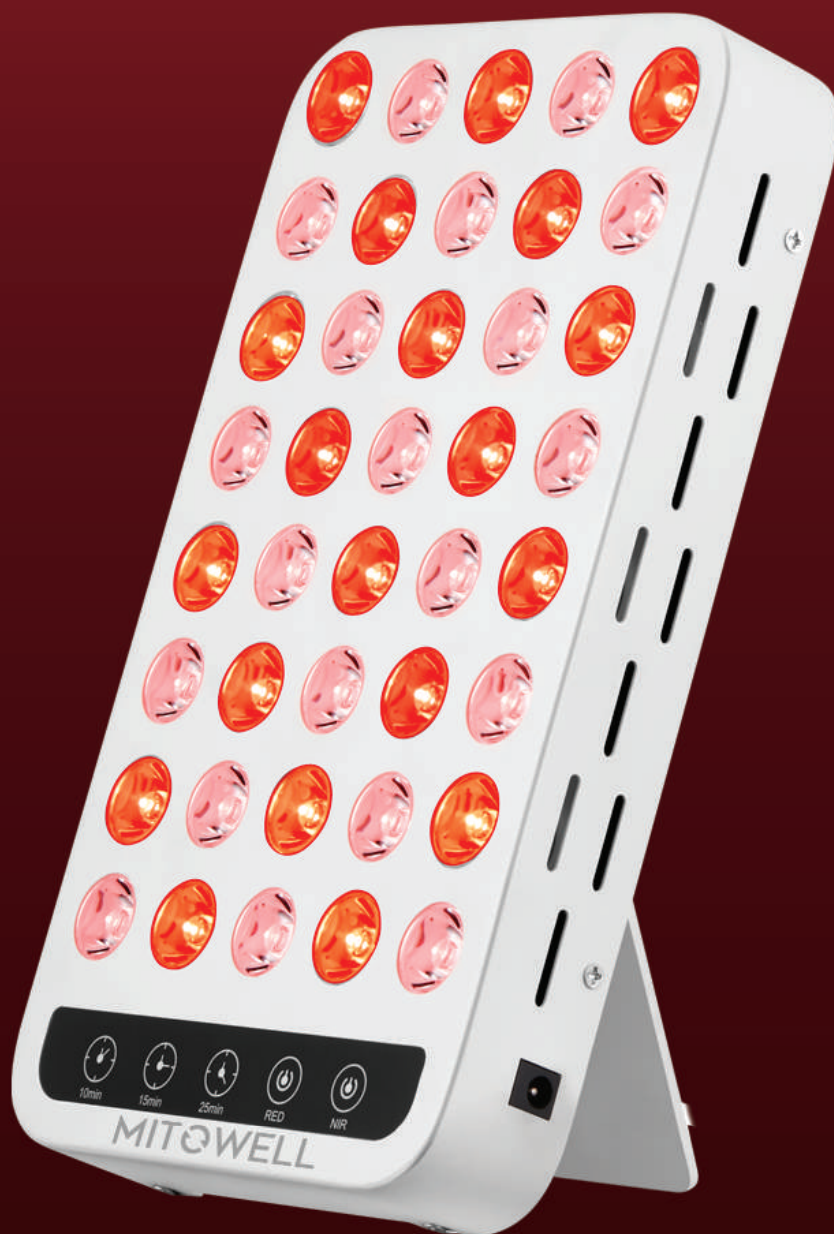
M200

# CONTENTS

INDICATIONS FOR USE .....	3
INSTRUCTIONS FOR USE .....	3
WHAT IS RED LIGHT THERAPY? .....	5
TECHNICAL SPECIFICATIONS .....	8
FEATURES AND CONTROLS .....	8
COMPOSITION OF THE SET .....	9
CONTRAINDICATIONS .....	9
POTENTIAL SIDE EFFECTS .....	9
SAFETY INFORMATION .....	10
DEVICE CLEANING AND CARE .....	11
DISPOSAL .....	11
WARRANTY / SERVICE .....	12
MANUFACTURER / IMPORTER / DECLARATION OF CONFORMITY .....	13
WARNINGS & DISCLAIMERS .....	13
END NOTES .....	14

# INTRODUCTION

## INDICATIONS FOR USE

Mitowell red light therapy (RLT) products are not intended to diagnose, treat or prevent any disease.

Our products are low-risk and intended for the wellness / fitness industry.

They are designed to support overall health by stimulating cellular processes. This product is not intended for illumination of indoor residential spaces.

If you are in doubt about the effects of the RLT on your health, you should always consult your doctor.

# INSTRUCTIONS FOR USE

## INSTALLATION & STARTUP

1. Remove the device from the packaging.
2. Remove all packaging material before using the device.
3. Before use, inspect the enclosure and power cord for visible damage.
4. Place the unit on a flat, level and stable surface.
5. Make sure that you connect the device to a mains voltage source that conforms to the specifications on the nameplate.
6. Connect the power cord in such a way as to prevent accidental tripping over it.


## EXPOSURE TIME & FREQUENCY

For optimum benefits, including effects on deep tissues, sessions of 5-25 minutes are recommended for each area you wish to treat with RLT, keeping the device 15-30 cm away from the area of treatment.

It is good practice to start RLT with shorter sessions and lengthen them depending on the body's response. A maximum of two to three areas can be exposed to RLT per day, with a recommended interval of 6 hours between sessions.


Time and distance may vary depending on individual response to RLT. Infrared light (IR light, wavelength 830nm + 850nm) affects deep tissues and red light (630nm + 660nm) is beneficial for the skin. There are many other benefits.

Both types of light can be used at the same time and this is recommended, as red and IR light complement each other.
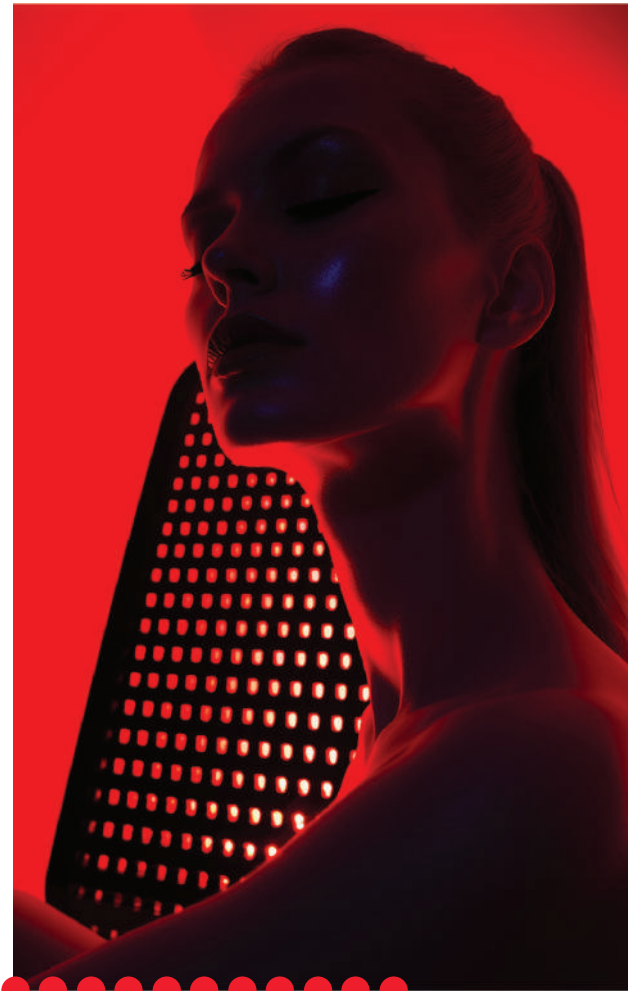
## EYE PROTECTION

Eye protection is not necessary, but staring directly into the LEDs is not recommended. If the light is uncomfortable for you, you can use the glasses included in the kit.
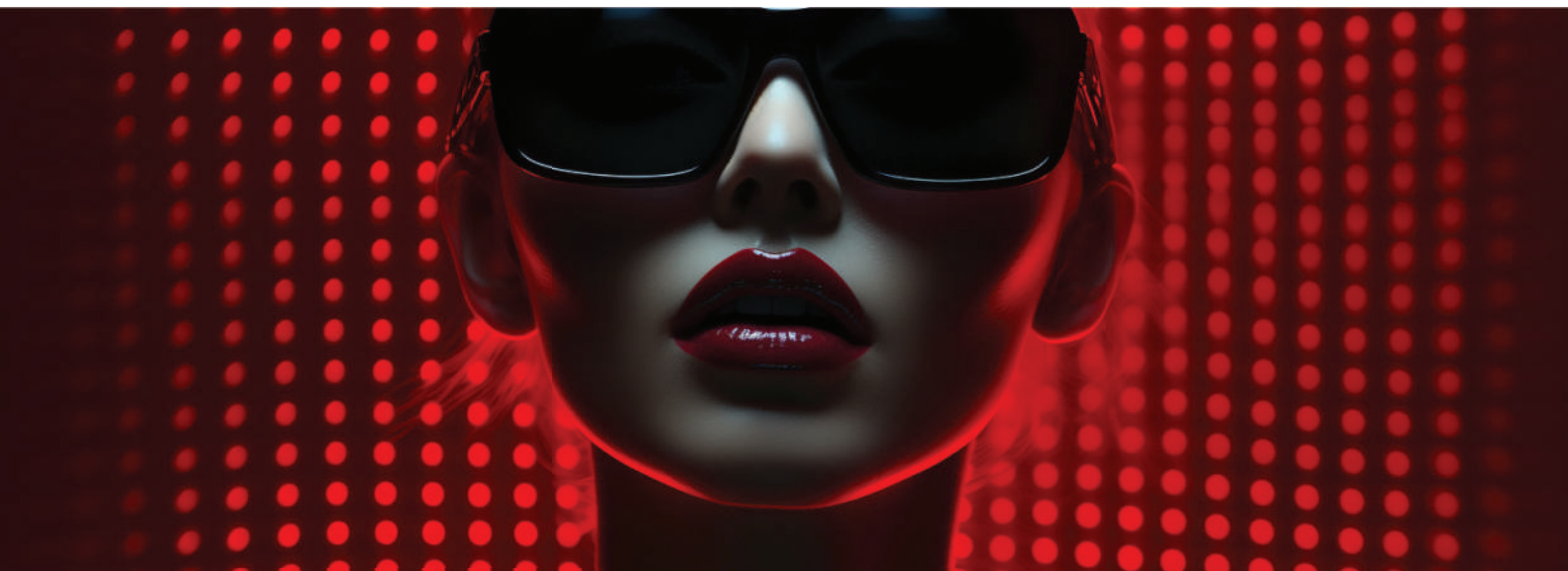
## GARMENTS

Make sure that the area of the body you are treating by RLT is completely exposed to the light and not covered by any fabric or clothing, or the benefits of RLT can be reduced.

# WHAT IS RED LIGHT THERAPY?

Light therapy is one of the oldest treatment methods known to mankind. The first uses of solar therapy were recorded in Egypt, and various forms of light therapy were also practised by the ancient Greeks, Chinese and Hindu.

There is no doubt that light affects the body biologically; in fact, a healthy body needs light.

Current research corroborates how different wavelengths of light affect the body at a cellular level.

All light wavelengths fall within a spectrum. Red and infrared (IR) light with a wavelength of 630 to 850 nanometres (nm) is particularly beneficial.

There are two key mechanisms responsible for the effectiveness of RLT:

### **Stimulation of ATP production**


RLT affects the mitochondria in cells, which are the 'power plants' that produce energy in the form of ATP. With red light irradiation, the production of ATP is increased, resulting in better performance of cells and the entire body.

### **Hormesis**

Short-term red light irradiation of the skin can induce a low stress level in the cells, resulting in a reinforcement of their natural defence mechanisms. This in turn contributes to better immunity and overall health.




In practical terms, this means that RLT can improve cellular processes and boost immunity. It is a method that supports health at a cellular level. For more details, refer to our RLT Guide included with the RLT lamp.

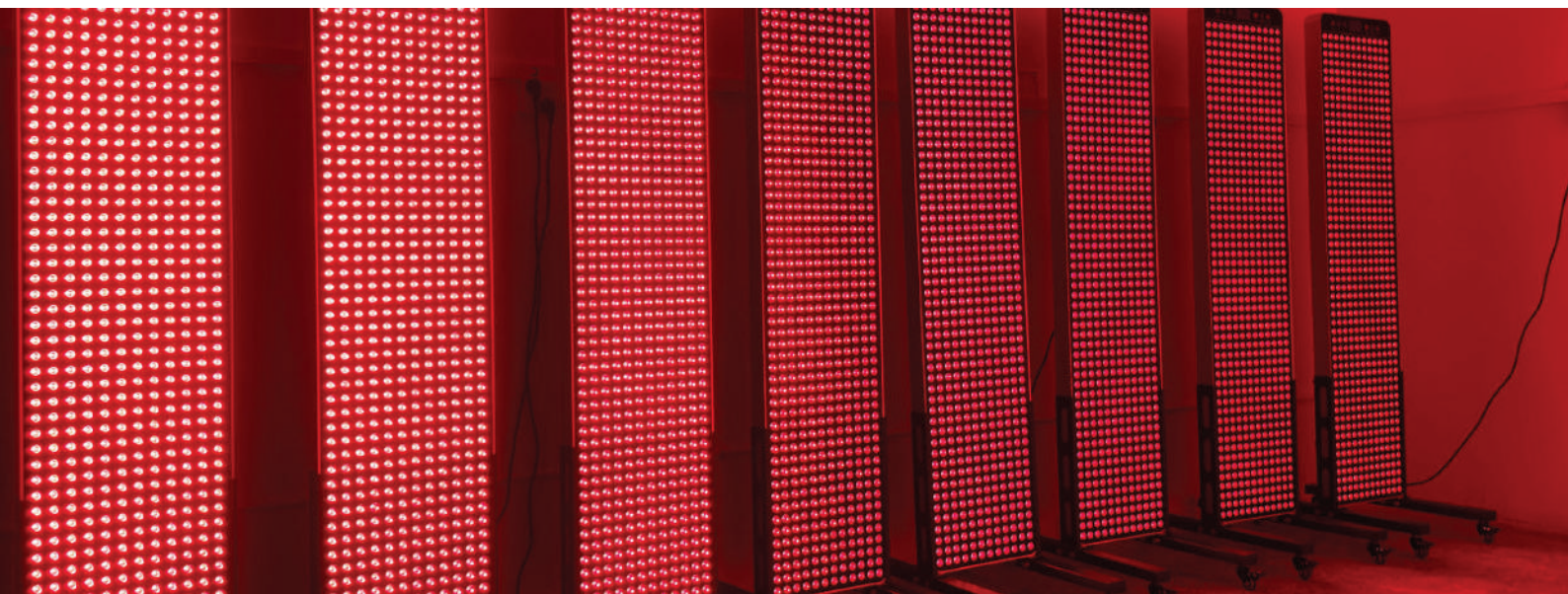
# TECHNICAL SPECIFICATIONS

Model	M200
Maximum LED output	200W
Power consumption	60W
Number of LEDs	40 pcs. DUAL-CHIP 5W
Wavelength	630, 660, 830, 850 nm (1:1:1:1)
Lens beam angle	30°
Service life	~ 50 000 hrs.
Irradiance, 3/15/30 cm distance (SOLAR METER)	~186/92/69 mW/cm <sup>2</sup>
EMF emission at 15 cm	0.0 V/m and 0.0 μT
Smart functions	3 light modes
1	● RED: 630, 660 nm
2	● NIR: 830, 850 nm
3	● RED&NIR: 630, 660, 830, 850 nm
Operating voltage	DC: 24V/2A
Stand mounting	No
Flicker	No
Noise (dB) at 15 cm	<53
Fans	1
RLT lamp dimensions	250X140X35 MM
RLT lamp weight	1.2 kg
Packaging dimensions	31*26*11 CM
Weight in packaging	2 kg

# FEATURES AND CONTROLS
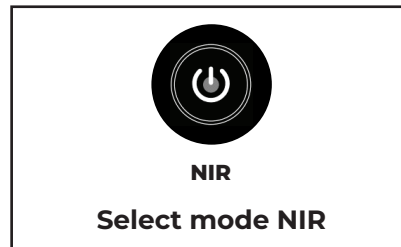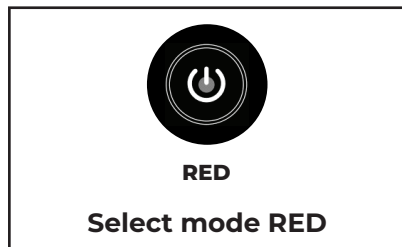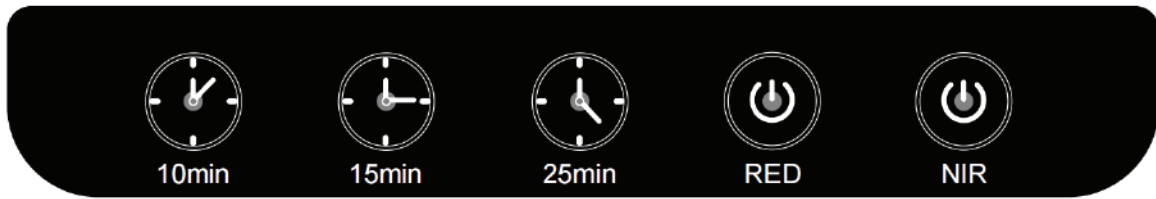
## COMPOSITION OF THE SET

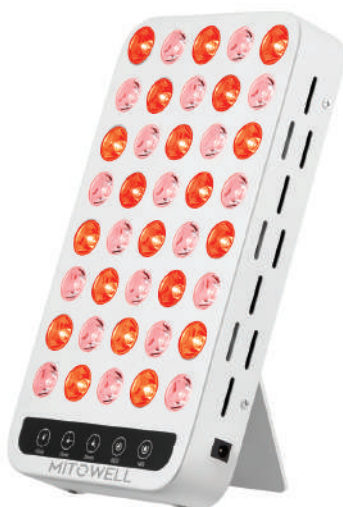Panel M200

Safety glasses

Power adapter

# CONTRAINDICATIONS

**The use of RLT is recognised as a safe and non-invasive method. Nevertheless, its use may not be recommended in certain circumstances.**

- Pregnancy and breastfeeding: it is recommended to avoid RLT for pregnant and breastfeeding women, due to the lack of conclusive scientific data.
- Photosensitisation and light sensitivity enhancing drugs: If you are using medications that increase light sensitivity, it is essential to consult your doctor before starting an RLT session.
- Eye disorders: People with ocular complaints or increased sensitivity to light should consult an ophthalmologist before attempting to undergo RLT.
- Cutaneous neoplasms: The use of RLT is not recommended for a diagnosis of skin cancer without prior consultation with an oncologist.
- Epilepsy: People with epilepsy must be cautious, as the therapy can induce an invisible flicker that may trigger a seizure.
- Skin burns: For fresh burns, RLT may not be recommended. Extra heat input can make the skin condition worse.
- Open wounds: RLT can irritate open wounds and affect the clotting process, which is not advisable.
- Fever: With elevated body temperature, the additional heat input generated by RLT may be inappropriate and aggravate the condition.
- Tattoos: People with tattoos, especially those with red ink, should be cautious as RLT may affect the colouring of the tattoo.

# POTENTIAL SIDE EFFECTS

**RLT is usually well tolerated and minimally invasive, but like any form of therapy, it can have some side effects:**

- Temporary redness at the irradiated site is one of the most common side effects.
- Dryness and irritation: After RLT, some people may experience dryness or irritation of the skin.
- Eye discomfort: Although red light is generally considered safe for the eyes, prolonged exposure can cause slight discomfort.
- Transient effects: On rare occasions, temporary symptoms such as a headache or a feeling of nausea may occur.

# SAFETY INFORMATION

To preserve the life and performance of the RLT LED panel, always store it in a cool and dry place. RLT devices shall only be operated in compliance with guidelines discussed in this Operating Manual. Carefully read and keep this Operating Manual for future reference, in a location accessible to other users.

## To ensure safety, do not do any of the following:

- Do NOT immerse the device in water
- DO NOT clean electrical components with water or other liquids
- DO NOT touch the device, plugs or wiring with wet hands
- DO NOT pull, carry or lift the device by its power cord
- DO NOT use the device if it has been dropped, left outdoors or has been exposed to water or other liquids.
- DO NOT use the device if it is damaged
- DO NOT use the device near other radiant heaters
- Do NOT modify or attempt to repair the device yourself
- DO NOT use or store the device at high temperatures that are or may be above 54°C (130°F)
- DO NOT use or store the device at a relative humidity above 93%
- DO NOT allow children to use the device without adult supervision
- DO NOT exceed 20 minutes of RLT exposure per 1 session unless otherwise advised by your doctor
- DO NOT cover the fans when using the device
- DO NOT connect the device to a source voltage other than that specified on the nameplate
- People with a limited physical, motor or mental capacity should NOT operate the device without the supervision of a competent person or without prior training in the safe use of the device.
- DO NOT leave the device connected to mains when not in use

# DEVICE CLEANING AND CARE


**We recommend cleaning the device as follows:**

- Unplug the device from the power supply and wait for it to cool down.
- Use only a damp microfibre cloth and make sure no water gets into the device
- Do not use strong cleaning products
- The device shall only be powered on again when it is completely dry.

---

## DISPOSAL

The crossed-out wheelie bin symbol shown on products or accompanying documentation indicates that end-of-life or defective electrical or electronic equipment shall not be disposed of with household waste.

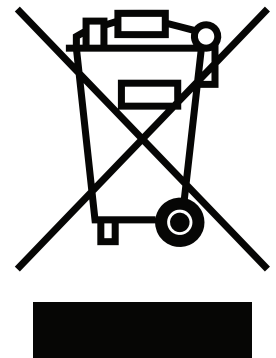The correct procedure when disposing of, reusing or recovering components is to take the device to a dedicated waste collection point, where it will be accepted free of charge.

In some countries, this product can be returned to the local distributor when purchasing another device. Proper disposal of the device can help to conserve valuable resources and avoid negative impact on health and the environment, which can be compromised by inappropriate handling of waste.

Contact your local authorities for details of your nearest waste collection point. Improper disposal of waste is subject to penalties under the relevant local regulations. Business users should contact their supplier for the terms and conditions of waste device return.

# WARRANTY / SERVICE

**Device designation:** RLT (red light therapy) LED lamp

**Model:** M200

**Lot number:** see the nameplate

The distributor of mitowell brand products, GWBL Sp. z o.o. (hereinafter referred to as the “Guarantor”) with its registered office at ul. Starowiejska 22B, 34-120 Andrychów, Poland provides a commercial warranty for this product under the following terms, conditions and applicability specified below:

- If the product is not in conformity with the relevant sales contract, the Buyer shall have legal remedies at the Seller’s side and expense, and these remedies shall not be affected by this commercial warranty.
- The Guarantor provides a 36-month warranty from the date of purchase of a new, unused product. The Guarantor warrants that the product offered for sale is free of defects.
- The Liability under this limited warranty covers defects arising from root causes inherent in the product as sold. Under this limited warranty, the Guarantor shall rectify physical defects free of charge. The rectification of defects or faults under this limited warranty will be done by repairing the product, or by replacing it or replacing the defective items with new, defect-free ones.
- Defects or damage to the product shall be reported within 7 days of emergence by e-mail to [warranty@mitowell.eu](mailto:warranty@mitowell.eu)
- The defects disclosed and reported during the limited warranty period will be rectified as soon as possible and within 21 days from the date of delivery of the defective product to the dealer or manufacturer’s service. The warranty repair period can be extended up to 30 days if the repair requires importing spare parts from abroad.

## Limited warranty exclusions:

- normal wear and tear related to routine operation of the included accessories that wear out in the course of normal use (e.g. eye protection glasses, electrical batteries, etc.); the warranty does not cover mechanical damage caused by misuse, dropping the device or other forms of physical damage.
- this limited warranty does not apply in the event of damage resulting from incorrect installation, maintenance or operation of the device.
- any unauthorised modification, repair or opening of the device will void this limited warranty.
- this limited warranty does not cover damage caused by mains overvoltage or connection to a mains with incorrect voltage.
- this limited warranty does not apply to damage resulting from extreme environmental conditions, such as extreme temperatures, humidity or exposure to chemicals.
- this limited warranty does not apply if the device has been operated in deviation from the accompanying Operating Manual.
- this limited warranty does not cover damage caused by natural disasters such as fire, flood, earthquake, etc.
- this limited warranty does not cover damage resulting from improper handling or transport of the device.
- this limited warranty does not apply to products that were purchased second-hand or identified as “outlet” products at the time of purchase.

# MANUFACTURER / IMPORTER / DECLARATION OF CONFORMITY

## Manufacturer:

Gbeing Smart Ag Co., Ltd.  
2nd Floor, Building 5, No. 226,  
Guangtian Road, Songgang, Tangxiayong,  
Bao'an District, Shenzhen, China

**Gbeing Smart Ag Co., Ltd.** hereby declares that the panel specified this manual comply with the essential requirements and other relevant provisions of the **EMC 2014/30/EU, RoHS 2015/863 and (EU)2017/2102 amending Annex II to 2011/65/EU.**

**Gbeing Smart Ag Co., Ltd.** hereby declares that the power adapter specified this manual comply with the essential requirements and other relevant provisions of the **EMC 2014/30/EU, LVD 2014/35/EU and RoHS 2011/65EU and its amendment EU 2015/863.**

## Declaration of conformity to be downloaded:

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## WARNINGS & DISCLAIMERS

Before you decide to begin RLT or any other form of treatment, we recommend consulting a qualified medical practitioner. Some people may be more sensitive to light therapy and adverse side effects may occur. The manufacturer does not guarantee the efficacy of RLT nor is it liable for its possible side effects.

Use RLT devices at your own risk. All health and treatment decisions should be made in consultation with a doctor or other qualified medical professional. The manufacturer who has issued this publication is not authorised to provide medical advice and the information contained in this manual is not intended for any diagnosis, treatment or prevention of any disease.

If you experience discomfort during irradiation with the LED panel or have doubts about RLT, stop using the device immediately and contact your doctor.

# MITOWELL

 [www.mitowell.eu](http://www.mitowell.eu)

 [info@mitowell.eu](mailto:info@mitowell.eu)

## USER MANUAL AND VIDEO TUTORIALS

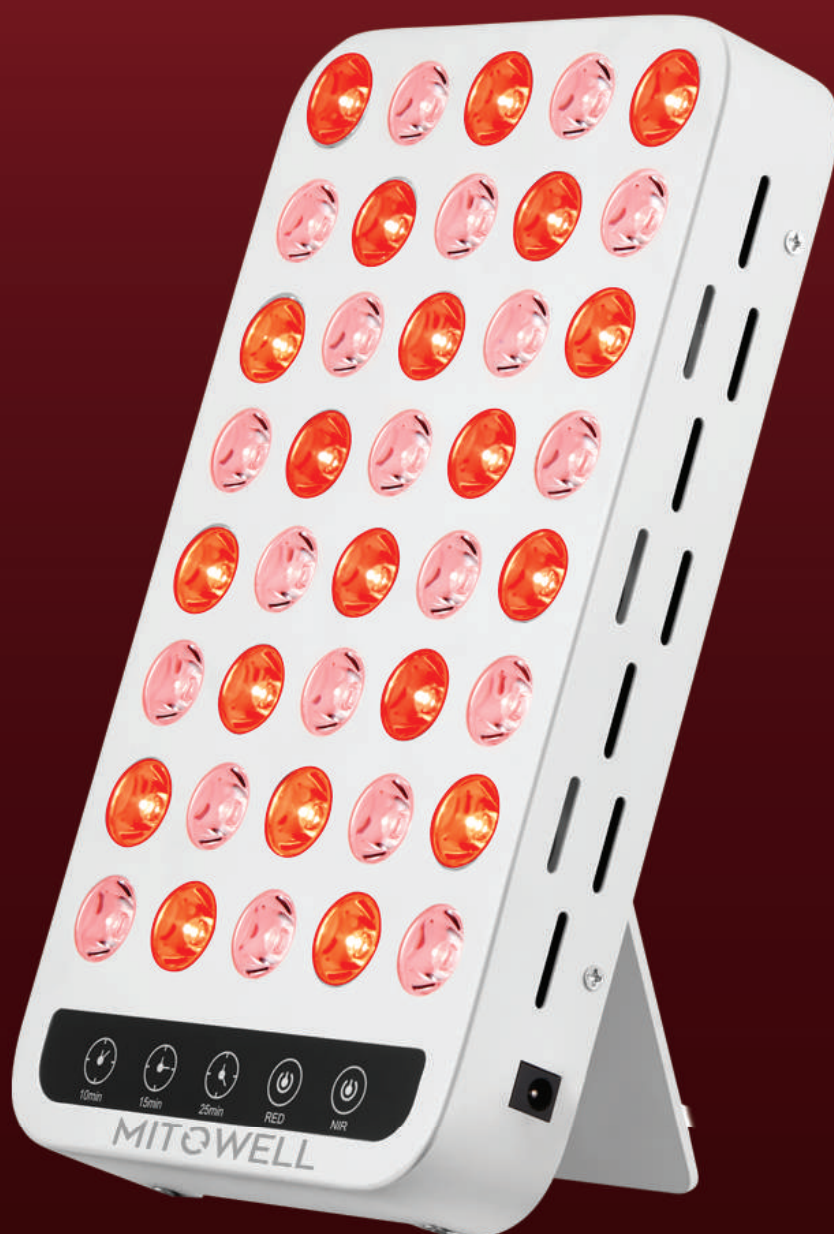
Find the user manual and video tutorials here:

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