
PL

EN

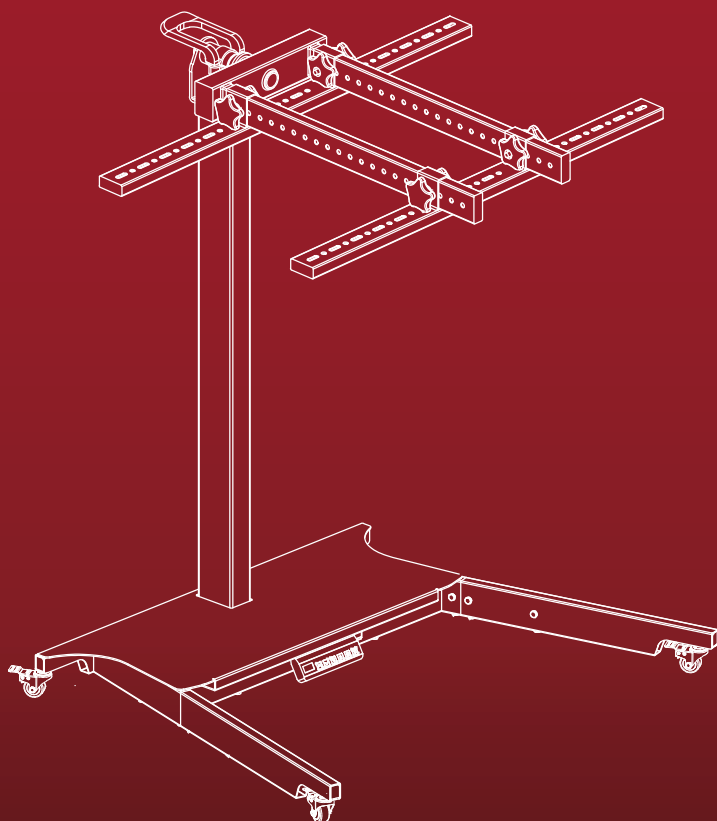MITOWELL

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ET01

EST-S

SPIS TREŚCI

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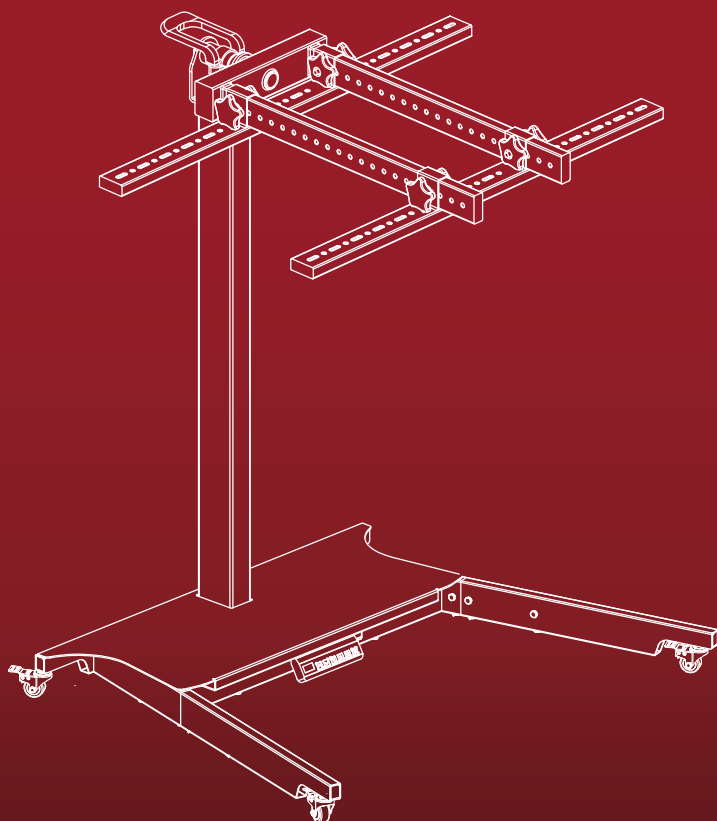
EST-S

TABLE OF CONTENTS

INTRODUCTION	3
OPERATING INSTRUCTIONS	3
INCLUDED ACCESSORIES EST-S	4
ASSEMBLY OF EST-S	5
INCLUDED ACCESSORIES ET01	7
ASSEMBLY OF ET01	8
MOUNTING THE PANEL ON THE STAND	10
ROTATION	11
REMOTE CONTROL	11
CONTROL PANEL	12
ERROR CODES AND SOLUTIONS	13
WARRANTY AND SERVICE	13
MANUFACTURER & IMPORTER	14
DECLARATION OF CONFORMITY	14

INTRODUCTION

PRODUCT DESCRIPTION

The most versatile stand for light therapy sessions – available in two sizes to suit different user needs.

Thanks to adjustable height and a rotating head, you can use therapy panels in almost any position: lying down, sitting, or standing.


It is compatible with Mitowell panels. It allows for both vertical and horizontal mounting, increasing the functionality of the devices during therapy.

OPERATING INSTRUCTIONS

NOTES






Please read the following warnings carefully before assembly:

1. Do not mount the lamp if the stand or equipment is damaged or missing components. If you need spare parts, contact us for assistance.
2. Do not use this product for purposes other than those intended and described in the user manual.
3. This stand contains moving parts. Use it with care.
4. The manufacturer is not responsible for any damage resulting from improper use, including unauthorized modifications, incorrect installation, or failure to follow these instructions.

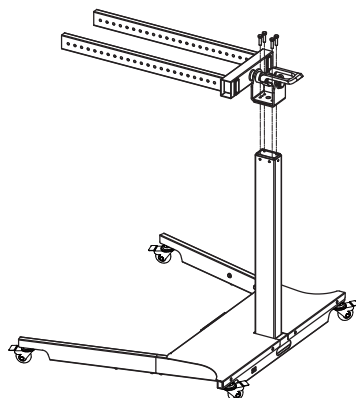INCLUDED ACCESSORIES EST-S

KIT CONTENTS

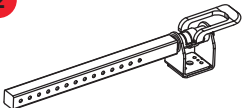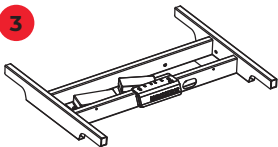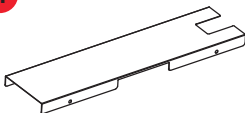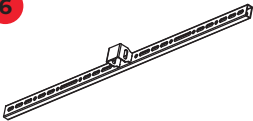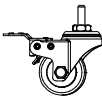1 	2 	3 
4 	5 	6 
7 	8 	9 
10 	11 	12 
13 	14 	15 

1 1 x Height-adjustable bracket

2 1 x Adjustable arm

3 1 x Control base

4 1 x Base cover plate

5 2 x Short mounting rails

6 2 x Long mounting rails

7 4 x Caster wheels

8 1 x Remote control

9 1 x Power cable

10 4 x Mounting screws

11 4 x Mounting nuts

12 1 x Assembly key


13 6 x Screw M6x12

14 12 x Screw M6x25


15 4 x Mounting bracket (angled)

ASSEMBLY OF EST-S


STEP-BY-STEP

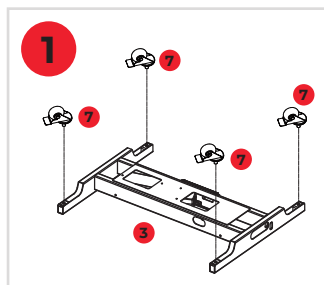
1. Attach the caster wheels (7) to the base (3).

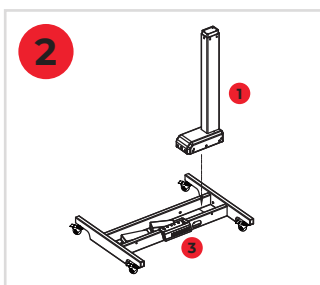
2. Insert the height-adjustable bracket (1) into the base (3).

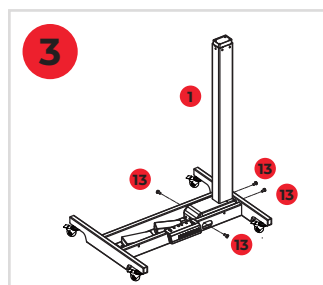
3. Secure the assembled structure using four screws (13) – see diagram.

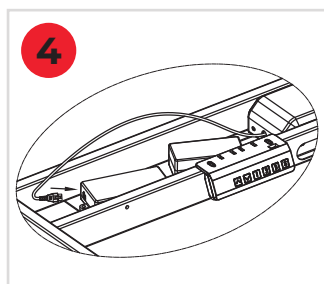4. Insert the cable from the bracket (1) into the base (3) – see diagram.

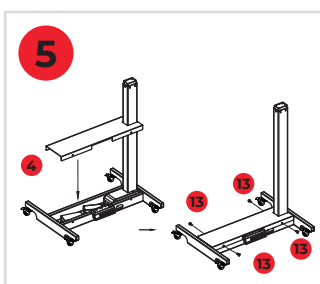
5. Place the cover plate (4) onto the base, then secure it with screws (13).

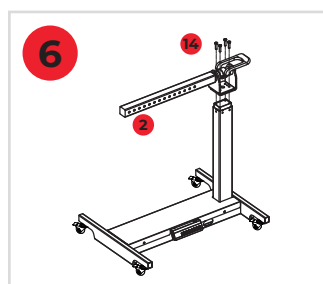
6. Attach the adjustable arm (2) to the bracket using screws (14).


PANEL INSTALLATION

MW600 PRO - 1 PCS

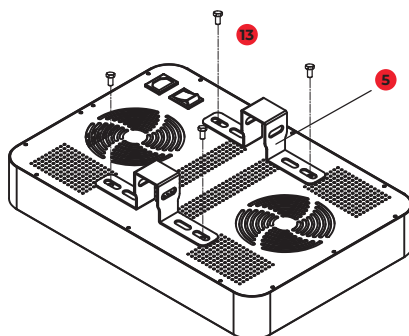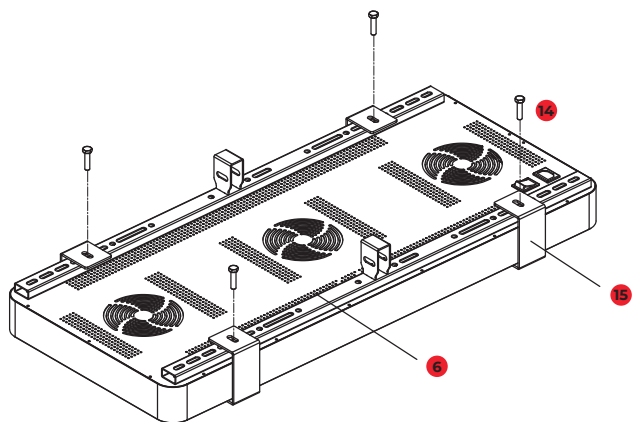
MW1500 PRO - 1 PCS

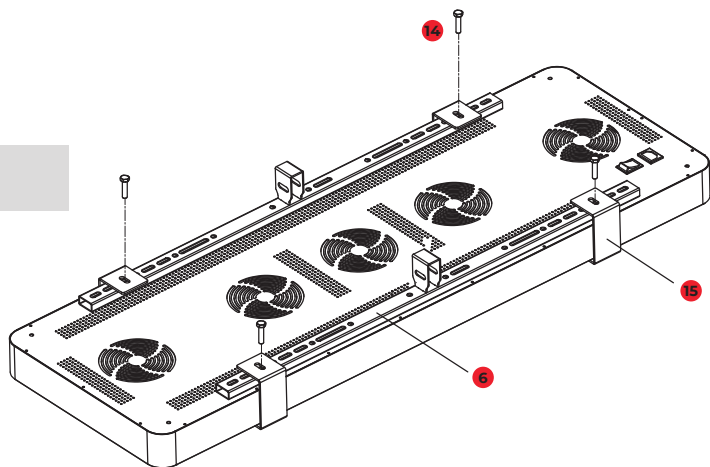
MW2000 PRO - 1 PCS

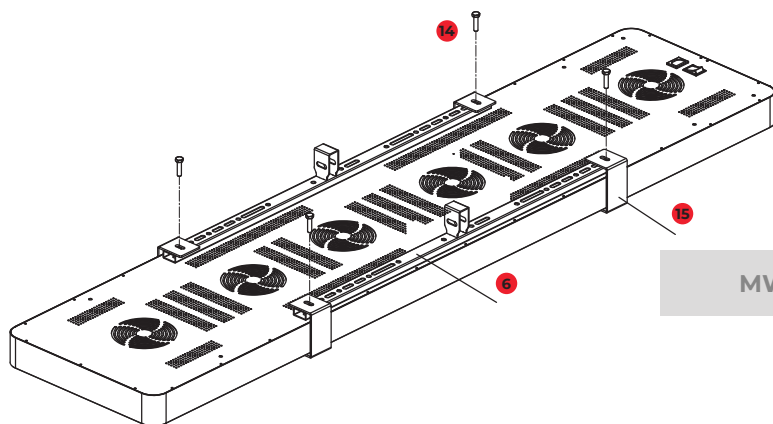MW3000 PRO - 1 PCS

INCLUDED ACCESSORIES ET01

KIT CONTENTS

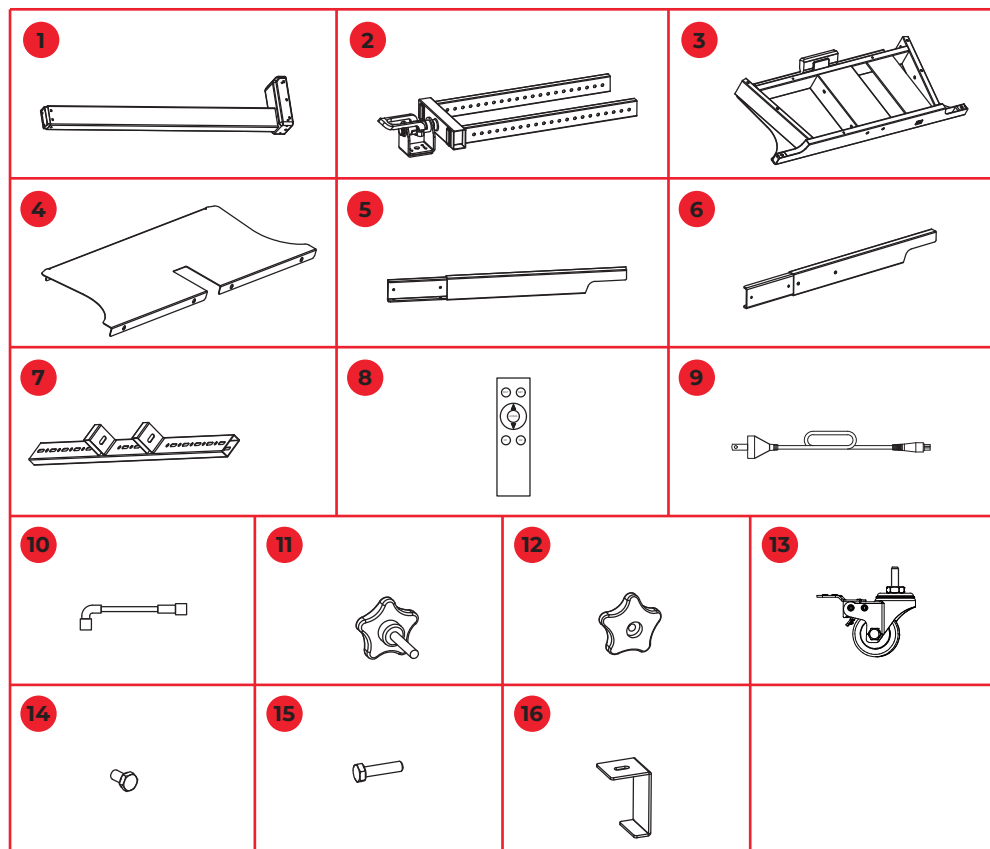1 1 x Height-adjustable bracket

2 1 x Adjustable arm

3 1 x Control base

4 1 x Base cover plate

5 1 x Right-side base mount

6 1 x Left-side base mount

7 2 x Mounting rails

8 1 x Remote control

9 1 x Power cable

10 1 x Assembly key

11 4 x Mounting screws

12 4 x Mounting nuts

13 4 x Caster wheels


14 16 x Screw M6x12

15 10 x Screw M6x25


16 4 x Mounting bracket (angled)

ASSEMBLY OF ET01


STEP-BY-STEP

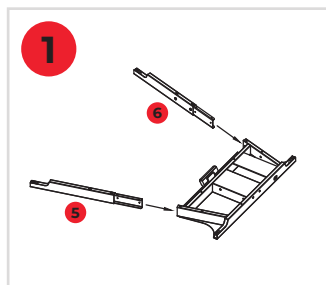
1. Insert the base mounts (5 and 6) into the base (3) – see illustration.

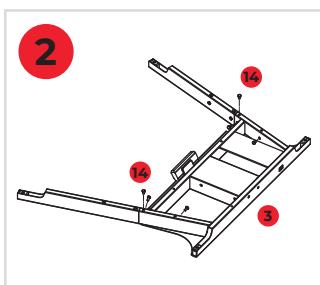
2. Secure the base mounts (5 and 6) using screws (14).

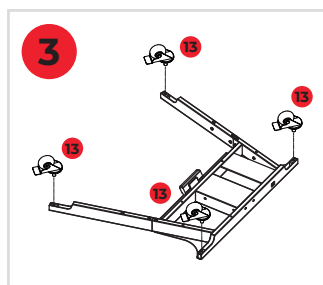
3. Attach the caster wheels (13) to the assembled base – see illustration.

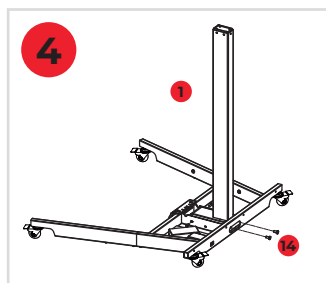
4. Insert the height-adjustable bracket (1) and fasten it to the base using screws (14).

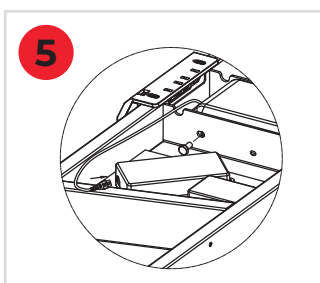
5. Insert the cable from the bracket (1) into the base (3) – see illustration.

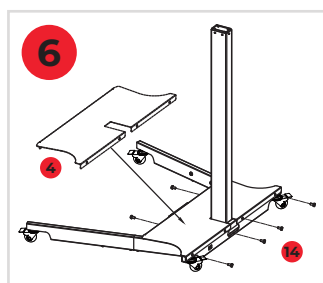6. Place the cover plate (4) onto the base and secure it using screws (14).


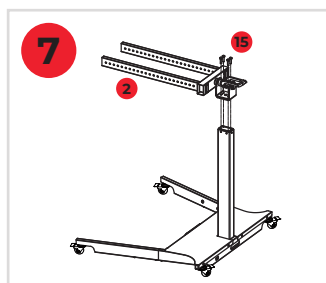
7. Attach the adjustable arm (2) and secure it to the bracket using screws (15).

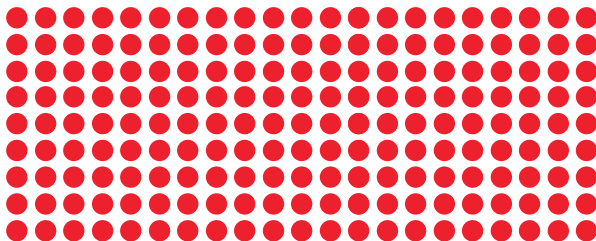
PANEL INSTALLATION

MW1500 PRO - 1 PCS

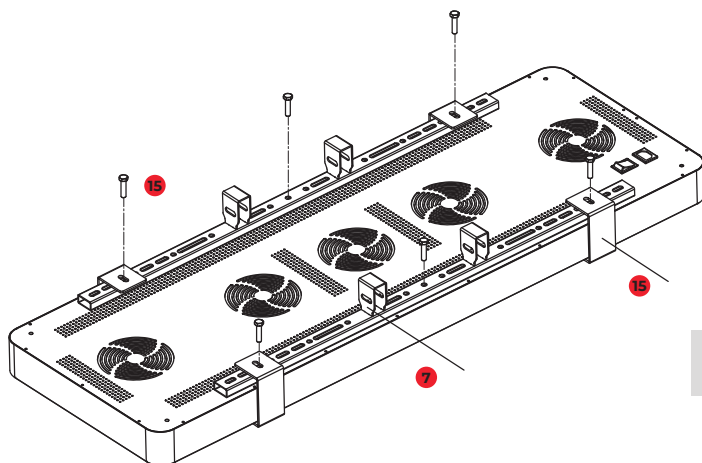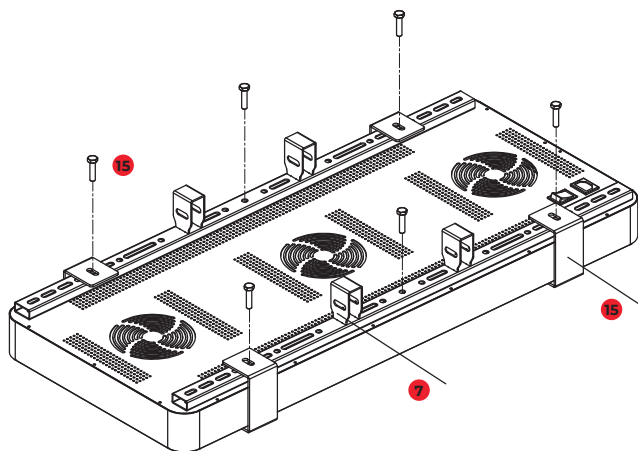MW2000 PRO - 1 PCS

MW3000 PRO - 1 PCS

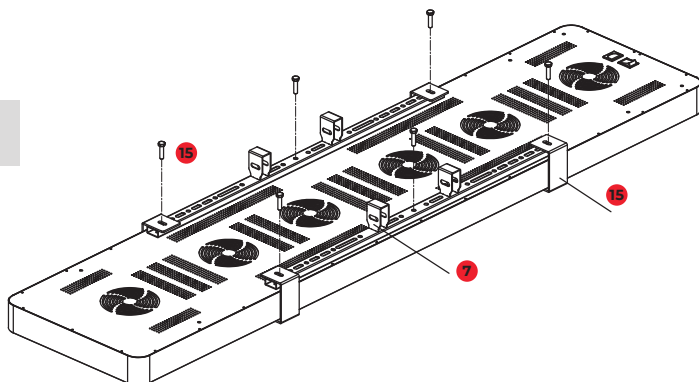MOUNTING THE PANEL ON THE STAND


MOUNTING PANELS ON THE EST-S STAND

EST-S is a solid and functional stand for Mitowell LED panels, designed to support up to 30 kg of weight.


It allows mounting of the following panels:

- MW600 Pro
- MW1500 Pro
- MW2000 Pro
- MW3000 Pro

It provides stability and flexibility for positioning the device during therapy sessions.

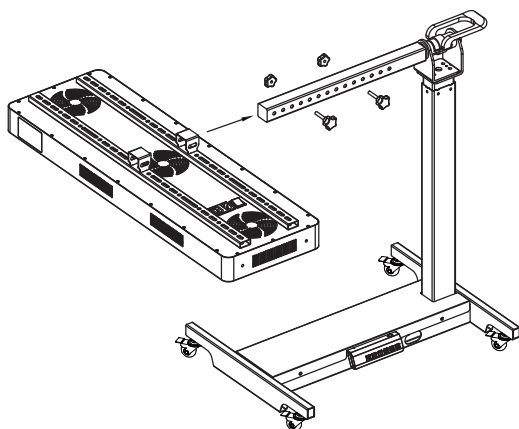MOUNTING PANELS ON THE ET01 STAND

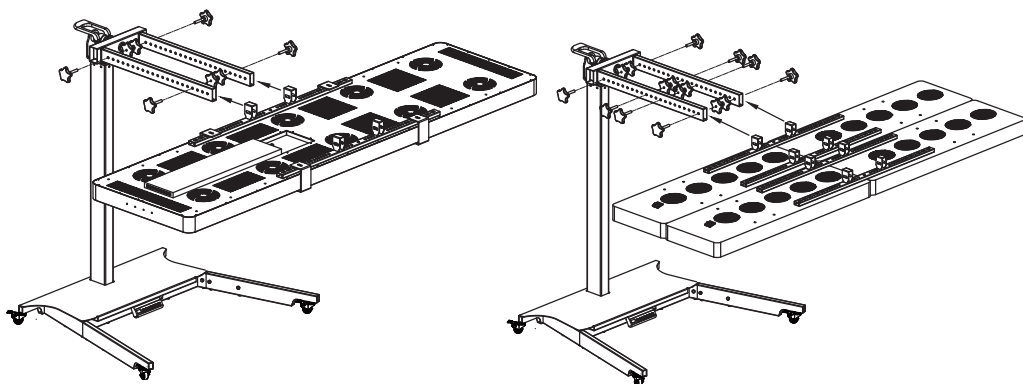
The ET01 is a durable stand with a maximum load capacity of 50 kg, enabling the mounting of one, two, or even four therapy panels. Thanks to its solid construction, it is ideal for both intensive home and professional use.

It allows mounting of the following panels:

- MW1500 Pro
- MW2000 Pro
- MW3000 Pro

Note: Brackets and screws required for mounting four panels are not included in the set – they must be purchased separately.

ROTATION

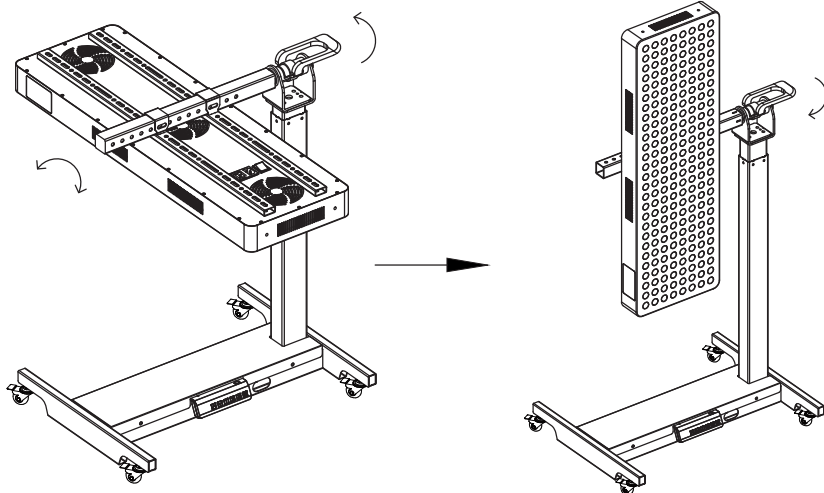1. Raise the stand to its highest position.
2. Hold the panel with one hand and lift the bracket upward.
3. Rotate the panel to your preferred position.
4. Press the bracket downward to lock it in place.

REMOTE CONTROL

Resetting

1. Press and hold the RST button for 3 seconds — you will hear a single beep indicating the reset process has started.
2. The stand will lower to the minimum height and bounce back up by 1 cm.
3. Two beeps will signal the reset is complete.

Remote Pairing


1. Unplug the stand's power cable.
2. Press and hold the SET button on the remote — do not release.
3. Plug in the power cable. When you hear the pairing signal, release the SET button.

Operation

- ▲ Press and hold to raise the stand; release to stop.
 - ▼ Press and hold to lower the stand; release to stop.
- HOME – Press to lower the stand to its minimum height. Press again to stop.
- M1 – Saved Height 1
- M2 – Saved Height 2


Saving Height

1. Move the stand to your desired height.
2. Press and hold M1 or M2 for 3 seconds — a beep will confirm the save.
3. To recall the saved height, press M1 or M2.

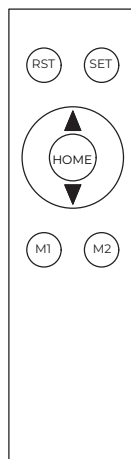CONTROL PANEL

PANEL OVERVIEW AND OPERATION

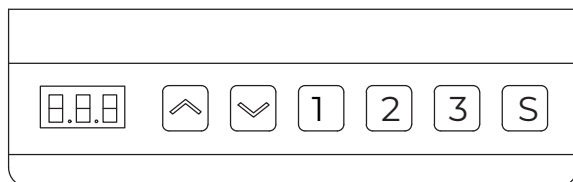First-time startup:

1. When the electric stand is powered on (100–240VAC), the buzzer emits a single beep and the controller displays the current height.
2. If the height displayed on the controller does not match the actual height, or if the stand was reassembled, a reset must be performed. Otherwise, the system can be used directly with factory settings.

Resetting

1. Press and hold the ▼ button; the stand will stop after 5 seconds.
2. “rst” will appear on the display — reset process begins.
3. Keep holding the ▼ button — the stand will slowly lower to its minimum height.
4. When it reaches the lowest point and bounces up by 1 cm, release the ▼ button to complete the reset.
5. In case of emergency, you can interrupt the reset anytime by pressing the ▼ button to stop the stand.

After an interrupted reset, the entire procedure must be repeated from the beginning.

Height adjustment

- Press and hold ▲ to raise the stand; release to stop.
- Press and hold ▼ to lower the stand; release to stop.

Saving height positions

1. Press the S button — “S–” appears on the display.
2. Press 1, 2, or 3 to save the current height under that button.
3. The control panel remembers up to three saved heights — even after disconnecting from power.
4. To recall a saved height, press the corresponding button — the stand will move to the saved position.
5. Press any button during movement to stop immediately.

Locking / Unlocking

1. Press and hold both ▲ and ▼ at the same time for 5 seconds.
2. “LOC” will appear on the display — the stand is now locked and cannot be moved.
3. To unlock, press and hold the same buttons again — the display will return to normal.

Energy saving mode

If no button is pressed for 1 minute, the stand enters energy-saving mode — the screen dims and power consumption is reduced. Press any button to exit this mode.

ERROR CODES AND SOLUTIONS

Error	Code	Solution
Overheat protection	HOT	Solution 1: No action required. The function will automatically recover after 18 minutes. Solution 2: Disconnect power, wait for the screen to turn off, then reconnect and resume use.
Reset not completed	rst	Press and hold the "RST" button on the control panel until the reset is fully completed.
Overload protection	E11	Remove heavy objects from the bracket. Wait for the alarm to stop before using the stand again.
Motor operation anomaly	E12	Turn off power, check the connection between the motor and controller. Turn the power back on.
Collision protection	E14	The stand will automatically reverse a specific distance. It should return to normal after the obstacle is cleared.
Motor lock alarm	E15	If there's a serious overload or the transmission mechanism is blocked, check whether any object is obstructing the movement.
Motor not connected	E16	Check if the motor is properly connected to the controller.

WARRANTY AND SERVICE

Product name: Panel Stand **Model:** ET01, EST-S

The distributor of Mitowell products, GWBL Sp. z o.o. (hereinafter referred to as the "Guarantor"), located at Starowiejska 22B, 34-120 Andrychów, Poland, provides warranty for this product under the following terms:

- In the event of non-compliance of the product with the sales agreement, the buyer has legal remedies at no cost. This warranty does not affect those rights.
- The Guarantor provides a 36-month warranty starting from the date of purchase of a new, unused product. The warranty covers material and manufacturing defects.
- The warranty covers defects arising from causes inherent in the product and does not include physical damage. Any defect outside the scope of this warranty will be treated as wear, misuse, or accidental damage.
- Any faults or damage must be reported within 7 days from the date they occur by email to: warranty@mitowell.eu
- Once a claim is submitted, the Guarantor will aim to process it as quickly as possible, within 21 days of receiving the defective product. This period may be extended up to 30 days if spare parts need to be imported from abroad.

The warranty does not cover:

- Normal wear and tear due to ongoing use of provided accessories during proper use.
- Mechanical damage caused by misuse, dropping the device, or other forms of physical damage.
- Damage resulting from improper installation, maintenance, or servicing of the device.
- Any unauthorized modifications — these void the warranty.
- Damage due to extreme environmental conditions (temperature, humidity, chemical exposure).
- Use of the device in a way not in accordance with the user manual.
- Damage caused by natural disasters (fire, flooding, earthquakes, etc.).
- Damage from improper transportation or handling.
- Products purchased as used or marked as "outlet" items.

MANUFACTURER & IMPORTER

Manufacturer:

Gbeing Smart Ag Co., Ltd.
2nd Floor, Building 5, No. 226,
Guangtian Road, Songgang, Tangxiayong,
Bao'an District, Shenzhen, China

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DECLARATION OF CONFORMITY

Gbeing Smart Ag Co., Ltd. hereby declares that the model described in this manual complies with the essential requirements and relevant provisions of applicable directives.

Declaration of conformity available at:

https://www.mitowell.eu/declaration/electric_stand.pdf

Made in China

WWW.MITOWELL.EU

INFO@MITOWELL.EU

USER MANUAL AND INSTRUCTIONAL VIDEOS

You can find the user manual and instructional videos here:

<https://mitowell.eu/user-manual>

SOCIAL MEDIA:

/MITOWELLEU

/MITOWELLEU

/@MITOWELLEU

/MITOWELLEU