
PL

EN
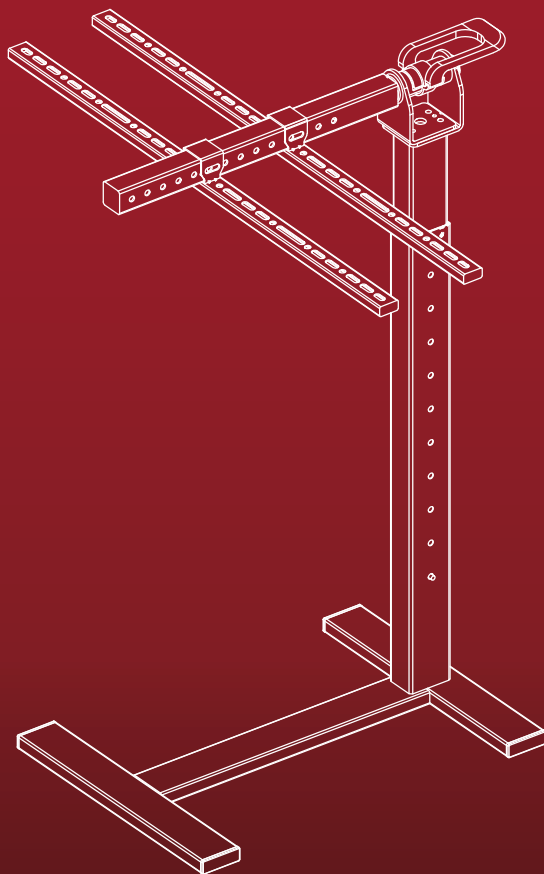
# MITOWELL

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MT02

# SPIS TREŚCI

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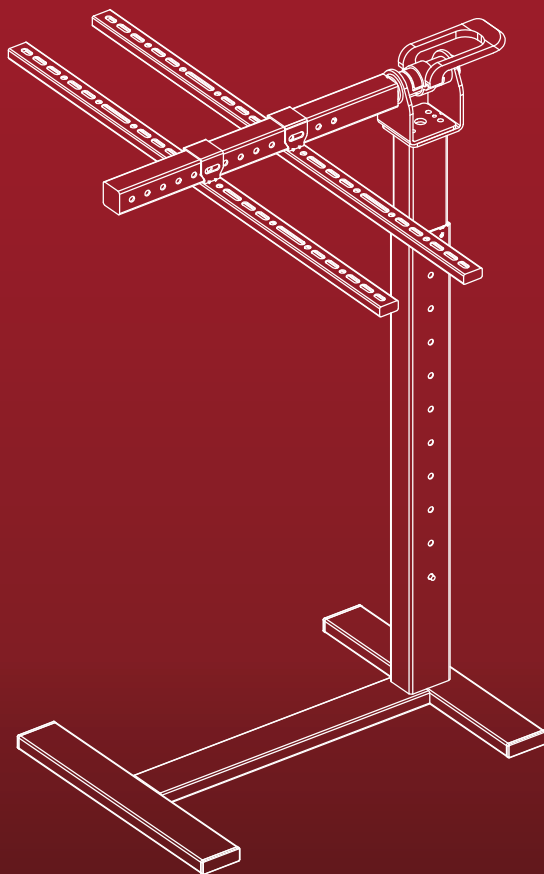
MT02

# TABLE OF CONTENTS

|                            |   |
|----------------------------|---|
| INTRODUCTION .....         | 3 |
| INSTRUCTIONS FOR USE ..... | 3 |
| INCLUDED ACCESSORIES ..... | 4 |
| ASSEMBLY .....             | 5 |
| WARRANTY AND SERVICE ..... | 9 |

# INTRODUCTION

## PRODUCT DESCRIPTION

A compact and functional stand allows for full 360° rotation, enabling precise positioning of the light during therapy.

Thanks to its lightweight and stable construction, it is perfect for both home and professional use.

The stand features smooth height adjustment from 945 mm to 1395 mm, allowing the panel's position to be adapted to individual needs.


# INSTRUCTIONS FOR USE

## NOTES


**Read the following warnings carefully before assembly.**

1. Mounting the lamp is not recommended if the stand or equipment is damaged, or if parts are missing from the set. If you need replacement parts, contact us for assistance.
2. The maximum load capacity of the stand is 16 kg.
3. Do not use this product for any purpose other than that described in the instruction manual.
4. The stand contains moving parts, so it should be used with caution.
5. The manufacturer is not responsible for any damage resulting from improper use of the product, i.e., unauthorized modifications, incorrect installation, or use not in accordance with its intended purpose.
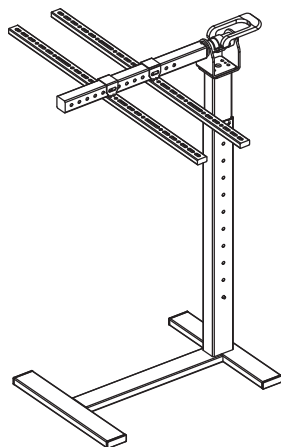
# INCLUDED ACCESSORIES

## PACKAGE CONTENTS

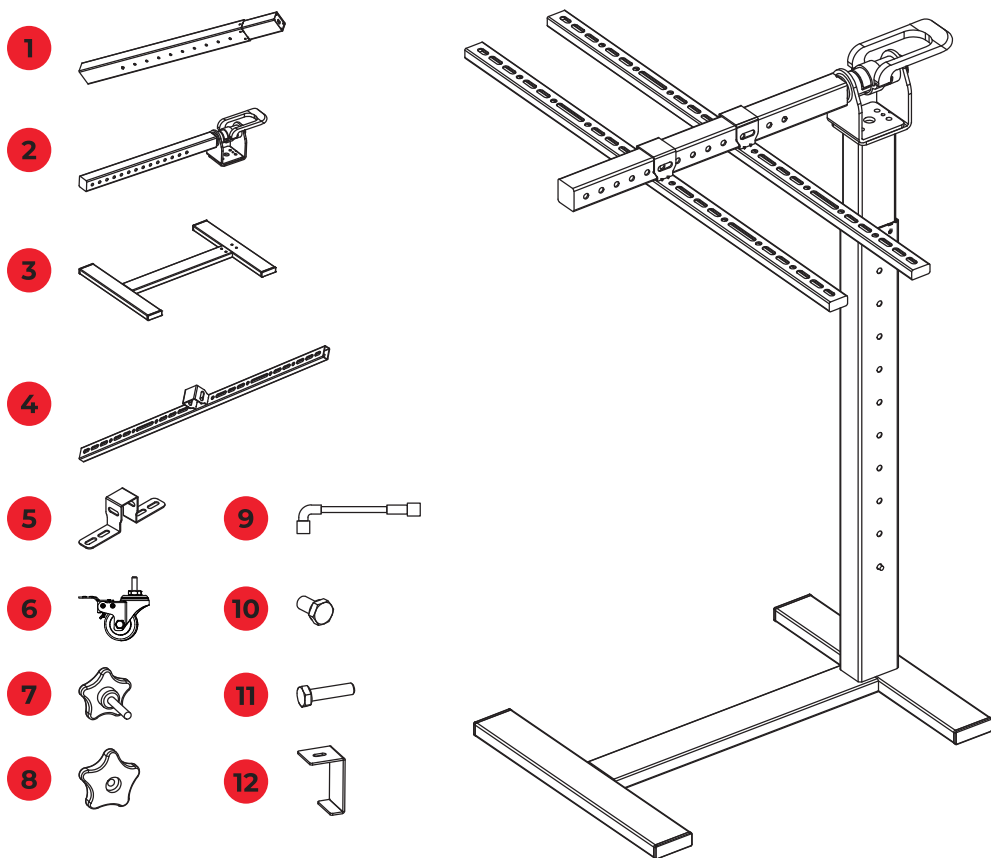
- |   |                       |    |                    |
|---|-----------------------|----|--------------------|
| 1 | Vertical support ×1   | 7  | Mounting screw ×3  |
| 2 | 360° Rotating arm ×1  | 8  | Mounting nut ×3    |
| 3 | Base ×1               | 9  | Assembly wrench ×1 |
| 4 | Adjustment bracket ×2 | 10 | M6x12 screw ×8     |
| 5 | Mounting bracket ×2   | 11 | M6x30 screw ×16    |
| 6 | Caster wheel ×4       | 12 | Corner bracket ×4  |


# ASSEMBLY

## Step by Step


**Note:** The base already has low-profile wheels installed, allowing basic mobility. If you want to increase mobility or height, install the additional wheels (6).  
If you do not want to install the larger wheels – skip to step 2.

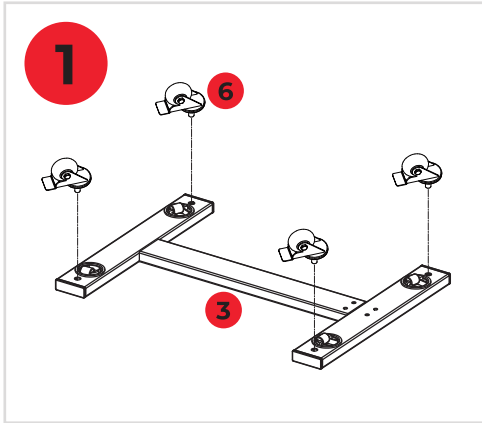
1. If you want to increase mobility or height, install the additional wheels (6). Check that they are properly tightened.

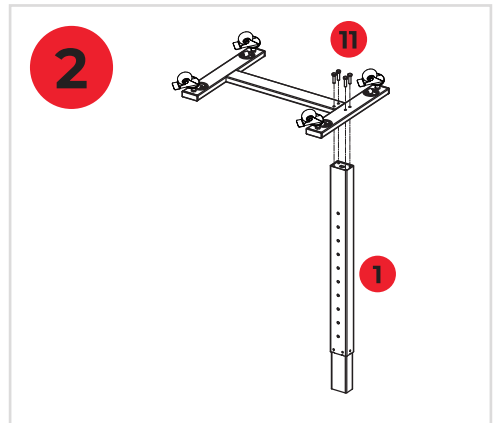2. Attach the vertical support (1) to the prepared base (3) using the screws (11).

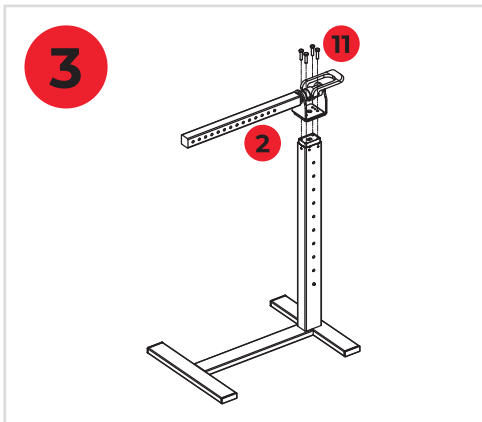
3. Attach the rotating arm (2) to the previously assembled structure using the screws (11) – see infographic.

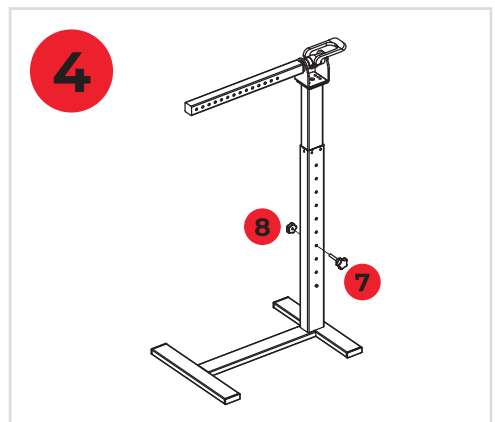4. Using the mounting screw (7) and mounting nut (8), set the stand according to your preferences. The stand allows height adjustment via locking in designated holes.


# ASSEMBLY


## PANEL CONFIGURATION

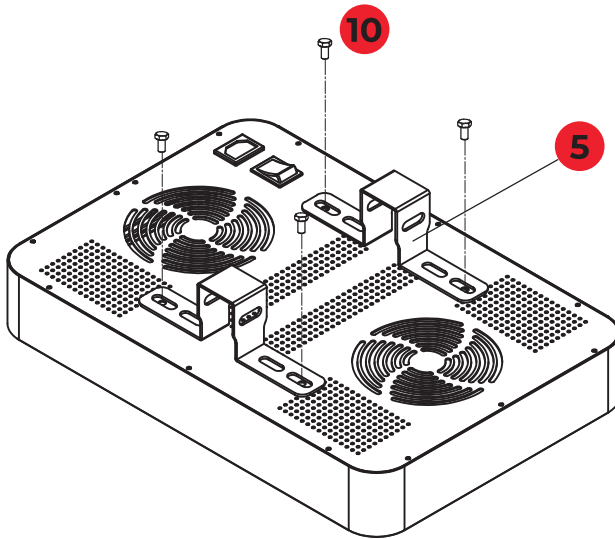
MW600 PRO

MW1500 PRO
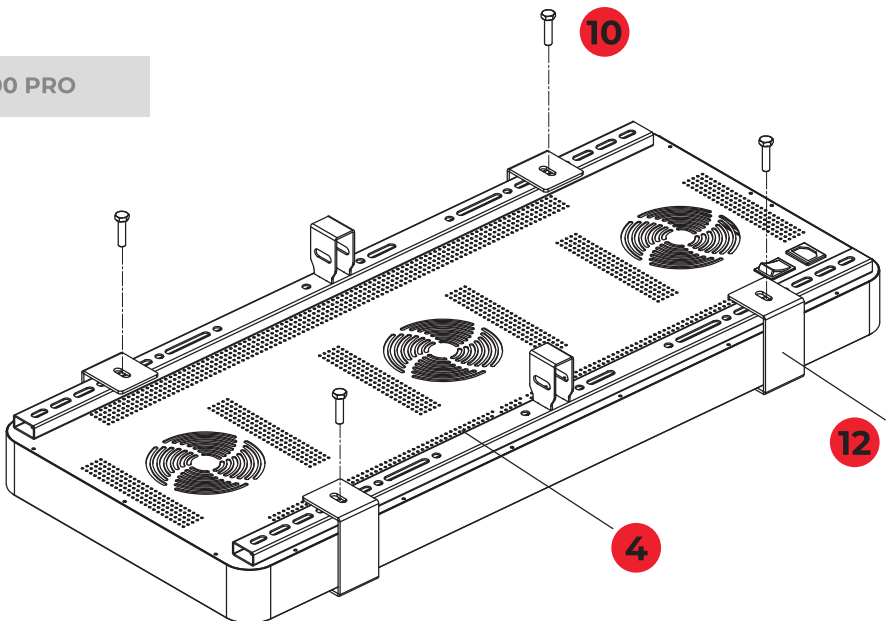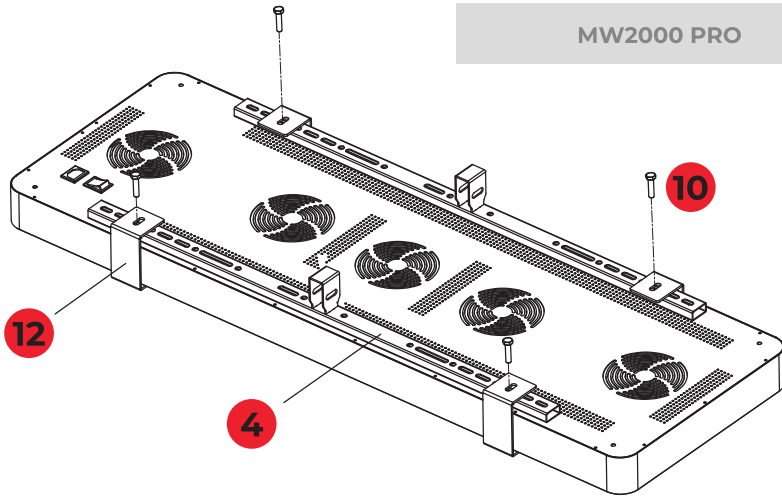
## VERTICAL / HORIZONTAL ADJUSTMENT OPTION
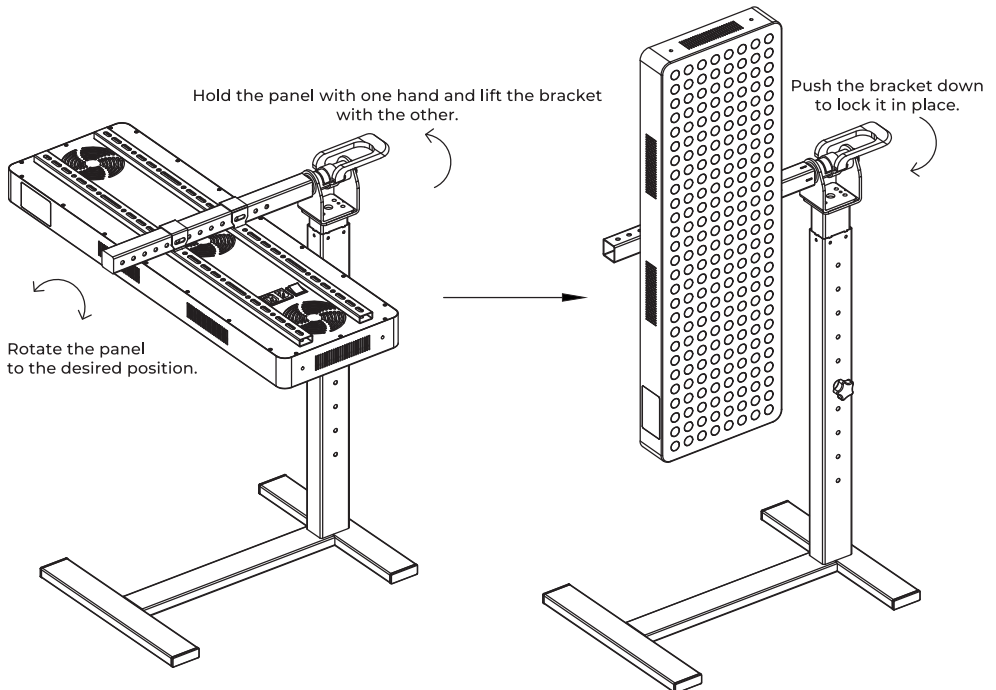
# WARRANTY & SERVICE

**Device designation:** Panel Stand

**Model:** MT02

The distributor of Mitowell brand products, GWBL Sp. z o.o. (hereinafter referred to as the "Guarantor") with its registered office at ul. Starowiejska 22B, 34-120 Andrychów, Poland provides a commercial warranty for this product under the following terms, conditions and applicability specified below:

- If the product is not in conformity with the relevant sales contract, the Buyer shall have legal remedies at the Seller's side and expense, and these remedies shall not be affected by this commercial warranty.
- The Guarantor provides a 36-month warranty from the date of purchase of a new, unused product. The Guarantor warrants that the product offered for sale is free of defects.
- The Liability under this limited warranty covers defects arising from root causes inherent in the product as sold. Under this limited warranty, the Guarantor shall rectify physical defects free of charge. The rectification of defects or faults under this limited warranty will be done by repairing the product, or by replacing it or replacing the defective items with new, defect-free ones.
- Defects or damage to the product shall be reported within 7 days of emergence by e-mail to [warranty@mitowell.eu](mailto:warranty@mitowell.eu)
- The defects disclosed and reported during the limited warranty period will be rectified as soon as possible and within 21 days from the date of delivery of the defective product to the dealer or manufacturer's service. The warranty repair period can be extended up to 30 days if the repair requires importing spare parts from abroad.

## Limited warranty exclusions:

- Normal wear and tear associated with the current use of supplied accessories that wear out during normal use are not covered by the warranty.
- The warranty does not apply to mechanical damage caused by improper use, device falls, or other forms of physical damage.
- The warranty does not apply in the case of damage resulting from improper installation, maintenance, or operation of the device.
- Any unauthorized modifications result in the loss of warranty.
- The warranty does not apply to damages resulting from extreme environmental conditions such as extreme temperatures, humidity, or exposure to chemicals.
- The warranty does not apply if the device was used in a manner not in accordance with the attached user manual.
- The warranty does not cover damages caused by natural disasters such as fires, floods, earthquakes, etc.
- The warranty does not apply to products purchased as used or labeled as "outlet".

**WWW.MITOWELL.EU**

**INFO@MITOWELL.EU**

## **USER MANUAL AND VIDEO TUTORIALS**

Find the user manual and video tutorials here:

**<https://mitowell.eu/user-manual>**

## **SOCIAL MEDIA:**

**/MITOWELLEU**

**/MITOWELLEU**

**/@MITOWELLEU**

**/MITOWELLEU**